

Trex Platform

ADVANCED

Veterinary Treatment & Exam Lift Table

This high-end treatment/exam table is versatile and is anchored by an adjustable powder-coated heavy-duty steel frame for maximum stability. Numerous configurations are available with manual or electric-powered units.

Primary options include battery or cord, mobility castors, integrated scale top, and steel versus synthetic tabletop surfaces.

FEATURES

- Height adjustment from 12-40 inches with manual or electric foot switch
- Stainless steel or synthetic top with marine edge and drain (51 x 23 in) and 14 tie-down hooks
- The optional synthetic top is made of a scratch-proof nonslip material. It is resistant to disinfectant and cleaning agents and eliminates the need for temporary mats
- Fixed or mobile platform with 4 optional braking castors
- Optional fixed or tiltable top (up to 11°) with one-hand operation supported by gas spring
- Optional integrated digital scale (0-220 lbs) with specialized animal weighing function. Backlit with USA or international units
- Battery-operated (220 hours working hours)
- 250lb load capacity



POWER VERSION



OPTIONAL: INTEGRATED SCALE



OPTIONAL: MOBILE CASTORS

Trex Platform Specifications



TECHNICAL SPECIFICATIONS

Table top length: 1300 mm (= 51.2 inch)
Table top width: 600 mm (= 23.6 inch)
Height adjustment: from ca. 310 – 1010 mm (=12.2 – 39.8 inch)

Maximum load capacity (electric): 120 kg (= 265 lbs)
Maximum load capacity (hydraulic): 100 kg (= 220 lbs)

Net weight: ca. 57 - 72 kg* (= 126 – 159 lbs)
Gross weight incl. packaging: ca. 76 - 91 kg* (= 168 – 201 lbs)

* Weight depending on width, accessories and other options

Height adjustment (with integrated scale): from ca. 335 –1045 mm (= 13.2 –41.2 inch)


Packaging dimensions: 1350 x 650 x 470 mm
(= 53.1 x 25.6 x 18.5 inch)

1350 x 650 x 700 mm (with integrated scale)
(= 53.1 x 25.6 x 27.6 inch)

Voltage (electrical version): 230V~, 50/60Hz, max. 1.5A
(valid for EU, CH, UK, AUS, CN)
120V~, 60Hz, max. 2.7A (USA)

Full load time: 10% (max. 2 min/20 min in operation)

Electrical data (battery): 24V DC, 2.9Ah

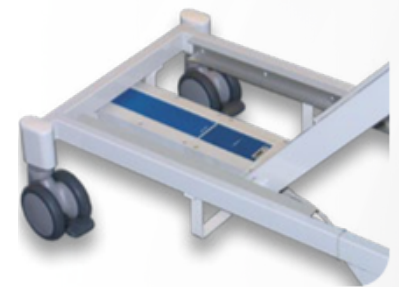


Battery power option



12 inch low height

Optional Scale Specifications:

- No. of Keys: 5
- 0-220 pounds in .1 pound units
- Auto Backlight: Yes, Selectable
- Auto Shutoff: Yes, Selectable
- LCD Height: 1 in
- USA or International measurement units

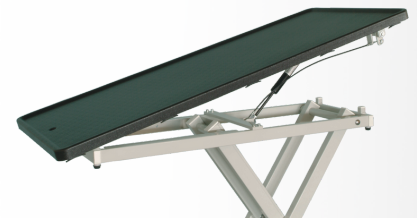

Gas supported tilt functions

OPTIONAL ACCESSORIES

- Quick-release tie-down cleat
- Infusion pole
- Adjustable side supports
- Narcotic bow
- Leg support with cleat
- Tilting stainless dish with grid
- Rail extension (200 mm, pair)
- Instrument tray (300x200 mm)
- Arm rest for ophthalmology
- Top extension (500x300 mm)
- Top extension (300x300 mm)
- Insert for OP-System



Optional Scale